

PAE212Hu71

Biotin-Linked Antibody to Elastase 3B (ELA3B)
Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: ELA3B, Human

Clonality: Polyclonal Conjugation: Biotin

Host: Rabbit


Immunoglobulin Type: IgG

Purification: Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

UOM: 100µg

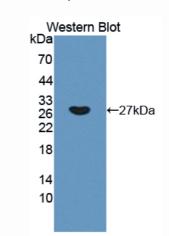

Sample: Recombinant ELA3B, Human

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant ELA3B (Val29~His270) expressed in E.coli.

Accession No.: RPE212Hu01

Sequence: The target protein is fused with N-terminal His-Tag and its sequence

is listed below.

MGHHHHHHSGS- VV NGEDAVPYSW PWQVSLQYEK SGSFYHTCGG SLIAPDWVVT AGHCISSSRT YQVVLGEYDR AVKEGPEQVI PINSGDLFVH PLWNRSCVAC GNDIALIKLS RSAQLGDAVQ LASLPPAGDI LPNETPCYIT GWGRLYTNGP LPDKLQEALL PVVDYEHCSR WNWWGSSVKK TMVCAGGDIR SGCNGDSGGP LNCPTEDGGW QVHGVTSFVS AFGCNTRRKP TVFTRVSAFI DWIFFTIASH



[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against ELA3B. It has been selected for its ability to recognize ELA3B in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:50-400

Immunocytochemistry in formalin fixed cells: 1:50-500

Immunohistochemistry in formalin fixed frozen section: 1:50-500

Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.