

RPH716Hu01 10 μ g

Recombinant Leucine Rich Repeat Kinase 2 (LRRK2)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Phe1479~Pro1729

Tags: N-terminal His-Tag

Accession: Q5S007

Host: *E. coli*

Subcellular Location: Cytoplasm. Membrane;
Peripheral membrane protein. Mitochondrion.

Purity: >95%

Endotoxin Level: <1.0EU per 1 μ g
(determined by the LAL method).


Formulation: Supplied as lyophilized form in PBS,
pH7.4, containing 1mM DTT, 5% trehalose, 0.01%
sarcosyl and preservative.

Predicted isoelectric point: 7.8

Predicted Molecular Mass: 31.0kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

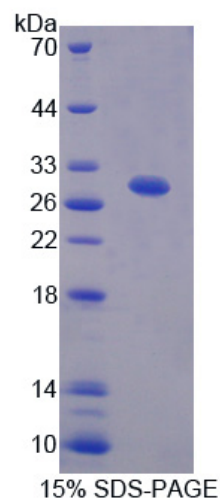

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

FP AIRDYHFVNA TEESDALAKL RKTIIINESLN FKIRDQLVVG QLIPDCYVEL EKIILSERKN
VPIEFVIDR KRLQLVREN QLQDENELP HAVHFLNESG VLLHFQDPAL QLSDLYFVEP
KWLCKIMAQI LTVKVEGCPK HPKGIISRRD VEKFLSKKRK FPKNYMSQYF KLEKFQIAL
PIGEEYLLVP SSLSDHRPVI ELPHCENSEI IIRLYEMPHYF PMGFWSRLIN RLEISPMYL
SGRERALRP

[REFERENCES]

1. Smith W.W., *et al.* (2005) Proc. Natl. Acad. Sci. U.S.A. 102:18676-18681.
2. West A.B., *et al.* (2005) Proc. Natl. Acad. Sci. U.S.A. 102:16842-16847.
3. Toff M., *et al.* (2005) Mech. Ageing Dev. 126:1201-1205.
4. Gloeckner C.J., *et al.* (2006) Hum. Mol. Genet. 15:223-232.