

**RPB503Hu02 100µg**

**Recombinant Dystrophin (DMD)**

**Organism Species: *Homo sapiens (Human)***

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Ile253~Lys597

**Tags:** N-terminal His Tag

**Subcellular Location:** Membrane, Plasma

**Purity:** > 90%

**Traits:** Freeze-dried powder

**Buffer formulation:** 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

**Original Concentration:** 400µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 42.9kDa

**Accurate Molecular Mass:** 45kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

IIGFRYRS LKHFNYDICQ SCFFSGRVAK GHKMHYPMVE YCTPTTSGED  
 VRDFAKVLKN KFRTKRYFAK HPRMGYLPVQ TVLEGDNMET PASSPQLSHD  
 DTHSRIEHYA SRLAEMENS N GSYLND SISP NESIDDEHLL IQHYCQSLNQ  
 DSPLSQPRSP AQILISLESE ERGELERILA DLEEENRN LQ AEYDRLKQQH  
 EHKGLSPLPS PPEMPTSPQ SPRDAELIAE AKLLRQHKGR LEARMQILED  
 HNKQLESQ LH RLRQLLEQPQ AEAKVNGTTV SSPSTSLQRS DSSQPMLLRV  
 VGSQTSDSMG EEDLLSPPQD TSTGLEEVME QLNN SFPSSR GRNTPGK

**[ IDENTIFICATION ]**

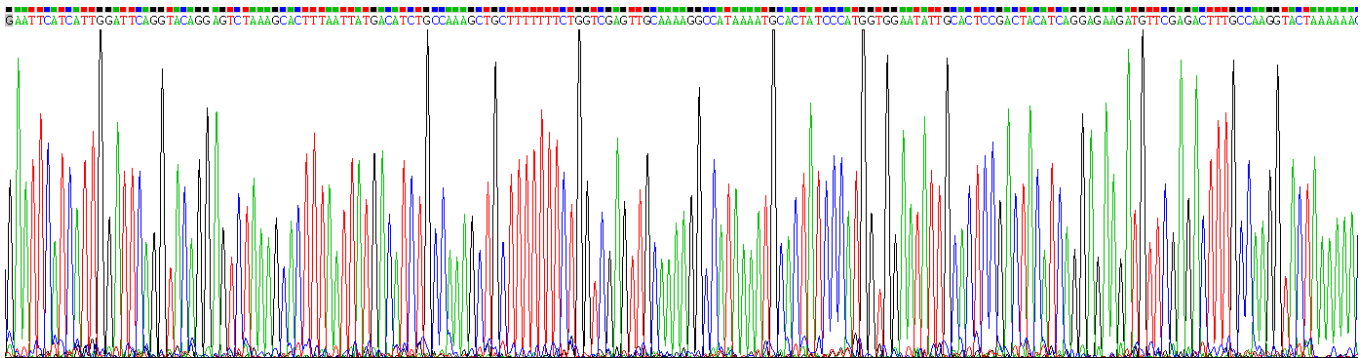

**Figure . Gene Sequencing (extract)**

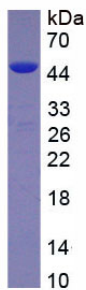


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.