

APA077Hu61 50μg

Active Interleukin 4 (IL4)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Eukaryotic expression.

Host: 293F cell

Residues: His25~Ser153 Tags: N-terminal His-tag

Purity: >95%

Endotoxin Level: <1.0EU per 1µg (determined by the LAL method).

Buffer Formulation: PBS, pH7.4, containing 5% Trehalose.

Original Concentration: 50µg/mL

Applications: Cell culture; Activity Assays; In vivo assays.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.3

Predicted Molecular Mass: 16.6kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.



[SEQUENCE]

HKCDIT LQEIIKTLNS LTEQKTLCTE
LTVTDIFAAS KNTTEKETFC RAATVLRQFY SHHEKDTRCL GATAQQFHRH
KQLIRFLKRL DRNLWGLAGL NSCPVKEANQ STLENFLERL KTIMREKYSK
CSS

[ACTIVITY]

The interleukin 4 (IL4,) is a cytokine that induces differentiation of naive helper T cells (Th0 cells) to Th2 cells. Upon activation by IL4, Th2 cells subsequently produce additional IL4 in a positive feedback loop. IL4 has many biological roles, including the stimulation of activated B-cell and T-cell proliferation, and the differentiation of B cells into plasma cells. It is a key regulator in humoral and adaptive immunity. IL4 induces B-cell class switching to IgE, and up-regulates MHC class II production. IL4 decreases the production of Th1 cells, macrophages, IFN-gamma, and dendritic cell IL12. The activity of IL4 is usually measured by a cell proliferation assay using TF-1 cells. TF-1 cells were seeded into triplicate wells of 96-well plates at a density of 20000 cells/well with 10% serum standard 1640 which contains various concentrations of recombinant human IL4. After incubated for 3 days, cells were observed by inverted microscope and cell proliferation was measured by Cell Counting Kit-8 (CCK-8). Briefly, 10 µl of CCK-8 solution was added to each well of the plate, then the absorbance at 450 nm was measured using a microplate reader after incubating the plate for 2-4 hours at 37 °C. Proliferation of TF-1 cells after incubation with IL4 for 3 days observed by inverted microscope was shown in Figure 1. Cell viability was assessed by CCK-8 (Cell Counting Kit-8) assay after incubation with rhIL4 for 3 days. The result was shown in Figure 2. It was obvious that rhIL4 significantly increased cell viability of TF-1 cells. The ED50 is 0.03 ng/ml.

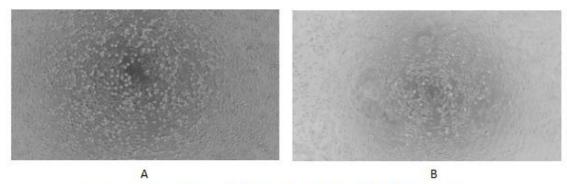


Figure 1. Cell proliferation of TF-1 cells after stimulated with IL4.

- (A) TF-1 cells cultured in 1640, stimulated with 0.1 ng/ml rhlL4 for 3 days;
- (B) Unstimulated TF-1 cells cultured in 1640 for 3 days.


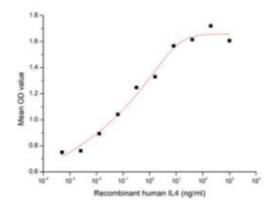


Figure 2. Cell proliferation of TF-1 cells after stimulated with rhIL4.

[IDENTIFICATION]

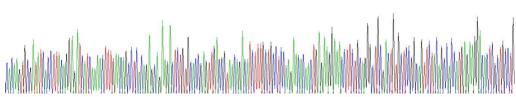

CALLEDGE MATTER CTRUCK TRANSCRIPT CREATER CONSIGNATION TO THE LIT VITOR TO LITERAL TRANSCRIPT TO THE ALS VITOR TO LITERAL TRANSCRIPT TO THE ALS VITOR TO THE CONTROL TO THE ALS VITOR OF ALS VITOR OF VEHICLE OF TREE LIT VITOR ALS VITOR OF VEHICLE OF TREE LIT VITOR ALS VITOR OF VEHICLE OF TREE LITERAL TO THE ALS VITOR OF VEHICLE OF THE CONTROL TRANSCRIPT TO THE ALS VITOR OF VEHICLE OF THE ALS VITOR OF THE ALS VITO

Figure 3. Gene Sequencing (extract)

Cloud-Clone Corp.

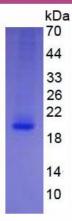


Figure 4. SDS-PAGE

Sample: Active recombinant IL4, Human

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.