

PAA281Hu02

Polyclonal Antibody to B-Cell Activation Factor Receptor (BAFFR)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 200μL

Cross Reactivity: N/A

Applications: IHC; ICC/IF

[IMMUNOGEN]

Immunogen: Recombinant BAFFR (Val96~Gln184 (Accession # Q96RJ3)) expressed in *E.coli*

Accession No.: RPA281Hu02

[APPLICATIONS]

Immunohistochemistry: 5-20μg/mL;

Immunofluorescence: 5-20μg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the

expiration date under appropriate storage condition.

[IDENTIFICATION]

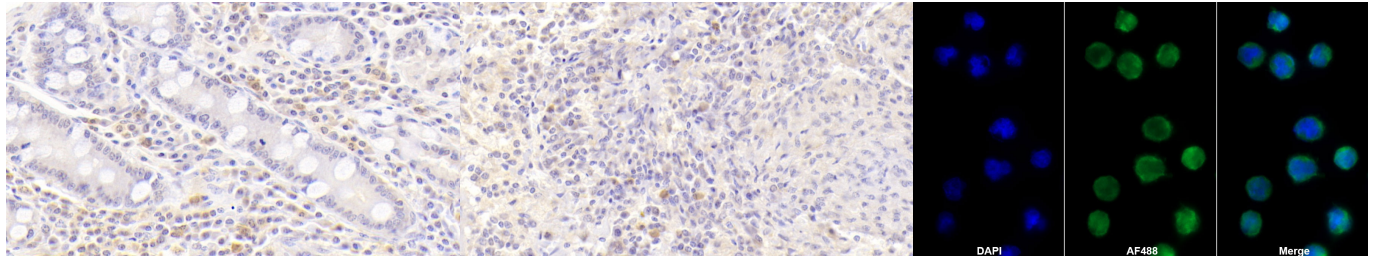

DAB staining on IHC-P;

Samples: Human Small intestine
Tissue;

Primary Ab: 20 μ g/ml Rabbit Anti-
Human BAFFR Antibody

Second Ab: 2 μ g/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody

(Catalog: SAA544Rb19)

DAB staining on IHC-P;

Samples: Human Spleen Tissue;
Primary Ab: 20 μ g/ml Rabbit Anti-
Human BAFFR Antibody

Second Ab: 2 μ g/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody

(Catalog: SAA544Rb19)

AF488 staining on IF;

Sample: Raji cell

Primary Ab: 20 μ g/ml Rabbit Anti-
Human BAFFR Antibody

Second Ab: 2 μ g/ml AF488-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody

(Catalog: SAA544Rb11)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.