

PAD071Ra01

Polyclonal Antibody to Carbonic Anhydrase IV (CA4)

Organism Species: *Rattus norvegicus* (Rat)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 100μL

Cross Reactivity: Human;Mouse;Porcine

Applications: WB; IHC; ICC/IF

[IMMUNOGEN]

Immunogen: Recombinant CA4 (Ser20~Ser220) expressed in *E.coli*

Accession No.: RPD071Ra01

[APPLICATIONS]

Western blotting: 0.01-2μg/mL;

Immunohistochemistry: 5-20μg/mL;

Immunocytochemistry: 5-30μg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no


obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

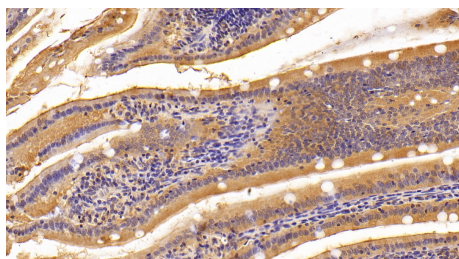


DAB staining on IHC-P; Sample: Rat


Small intestine Tissue; Primary Ab:
20ug/ml Rabbit Anti-Rat CA4 Antibody
Second Ab: 2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody (Catalog: SAA544Rb19)

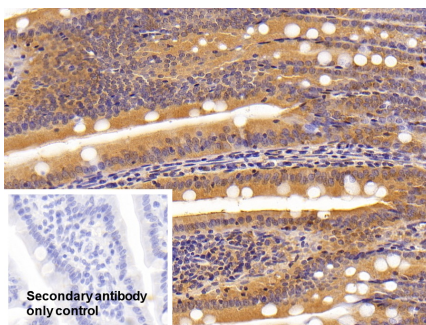

DAB staining on IHC-P; Sample: Rat
Small intestine Tissue Primary Ab:
20μg/ml Rabbit Anti-Rat CA4
Antibody Control: Used PBS instead of
primary antibody Second Ab: 2μg/ml
HRP-Linked Caprine Anti-Rabbit IgG
Polyclonal Antibody (Catalog:
SAA544Rb19)

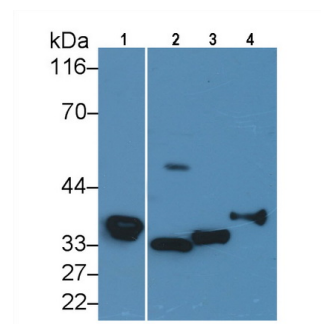

Western Blot; Sample: Lane1: Rat Lung
lysate; Lane2: Human Lung lysate;
Lane3: Mouse Lung lysate; Lane4:
Porcine Lung lysate
Primary Ab: 0.2μg/ml Rabbit Anti-Rat
CA4 Antibody
Second Ab: 0.2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

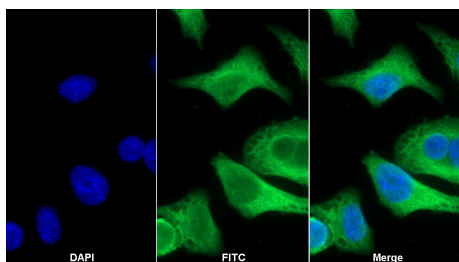

FITC staining on IF;

Sample: HeLa cell

Primary Ab: 30μg/ml Rabbit Anti-Rat
CA4 Antibody
Second Ab: 2μg/ml FITC-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb11)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.