

PAD195Hu01

Polyclonal Antibody to Collagen Type III Alpha 1 (COL3a1)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 1mL

Cross Reactivity: Mouse;Porcine

Applications: WB

[IMMUNOGEN]

Immunogen: Recombinant COL3a1 (Gly1059~Leu1466) expressed in *E.coli*

Accession No.: RPD195Hu01

[APPLICATIONS]

Western blotting: 0.01-2µg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]




Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

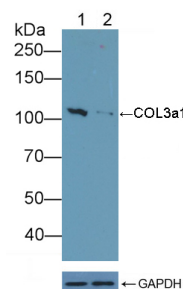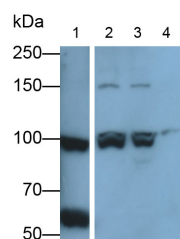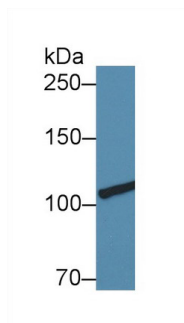

Western Blot; Sample: Porcine Spleen lysate; Western Blot; Samples: Lane1: Porcine

lysate;

Spleen lysate; Lane2: BXP3C cell

Knockout Varification:

Lane 1: Wild-type BXP3C cell lysate;

Primary Ab: 1µg/ml Rabbit Anti-Human

lysate; Lane3: HeLa cell lysate; Lane4:

Lane 2: COL3a1 knockout BXP3C cell

COL3a1 Antibody

RAW264.7 cell lysate;

lysate;

Second Ab: 0.2µg/mL HRP-Linked

Primary Ab: 0.3µg/ml Rabbit Anti-

Predicted MW: 111,138kd

Caprine Anti-Rabbit IgG Polyclonal

Human COL3a1 Antibody

Observed MW: 110kd

Antibody

Second Ab: 0.2µg/ml HRP-Linked

Primary Ab: 1µg/ml Rabbit Anti-Human

(Catalog: SAA544Rb19)

Caprine Anti-Rabbit IgG Polyclonal

COL3a1 Antibody

Antibody

Second Ab: 0.2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.