

PAB642Hu01

Polyclonal Antibody to Electron Transfer Flavoprotein Beta Polypeptide (ETFb)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 50µL

Cross Reactivity: Mouse;Rat;Porcine.

Applications: WB; IHC; ICC/IF

[IMMUNOGEN]

Immunogen: Recombinant ETFb (Ala2~Ile255) expressed in *E.coli*

Accession No.: RPB642Hu01

[APPLICATIONS]

Western blotting: 0.01-2µg/mL;

Immunohistochemistry: 5-20µg/mL;

Immunocytochemistry: 5-20µg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no

obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

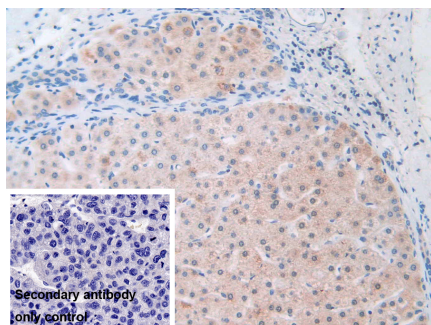

DAB staining on IHC-P; Sample:

Human Liver cancer Tissue Primary Ab: Samples: Human Breast cancer Tissue;


10µg/ml Rabbit Anti-Human ETFb

Antibody Control: Used PBS instead of primary antibody Second Ab: 2µg/ml

HRP-Linked Caprine Anti-Rabbit IgG

Polyclonal Antibody (Catalog:

SAA544Rb19)

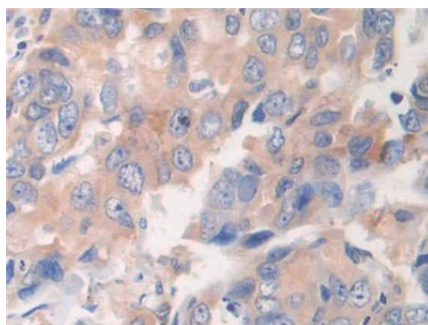

DAB staining on IHC-P;

Primary Ab: 10µg/ml Rabbit Anti-


Human ETFb Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

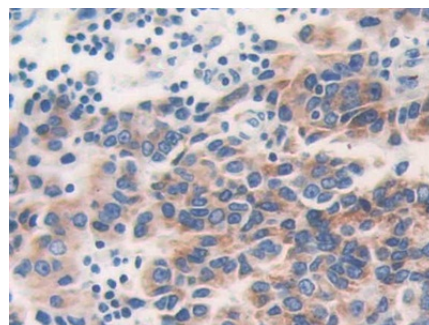

DAB staining on IHC-P;

Samples: Human Prostate cancer

Tissue;

Primary Ab: 10µg/ml Rabbit Anti-


Human ETFb Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

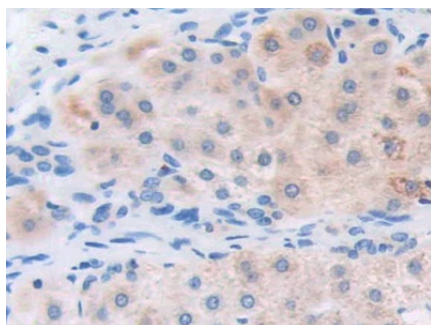

DAB staining on IHC-P;

Samples: Human Liver cancer Tissue;

Primary Ab: 10µg/ml Rabbit Anti-


Human ETFb Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

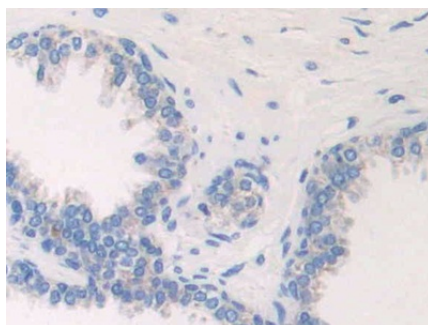

DAB staining on IHC-P;

Samples: Human Prostate Tissue;

Primary Ab: 10µg/ml Rabbit Anti-


Human ETFb Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

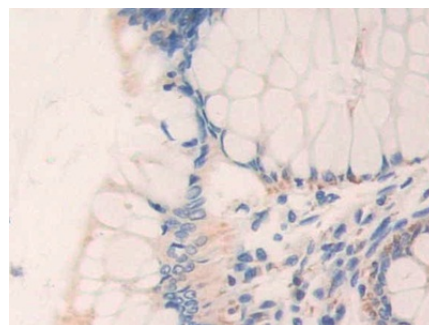

DAB staining on IHC-P;

Samples: Human Colorectal cancer

Tissue;

Primary Ab: 10µg/ml Rabbit Anti-



Human ETFb Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

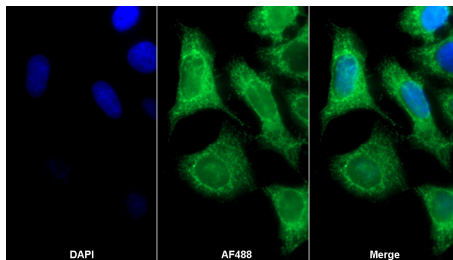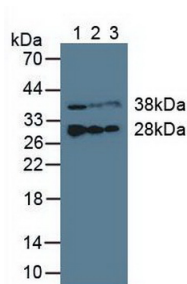

AF488 staining on IF;

Sample: HepG2 cell

Primary Ab: 20µg/ml Rabbit Anti-

Human ETFb Antibody

Second Ab: 2?g/ml AF488-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb11)

Western Blot; Samples: Lane1: Porcine

Liver lysate; Lane2: Porcine Heart

lysate; Lane3: Porcine Skeletal muscle

lysate;

Primary Ab: 3µg/ml Rabbit Anti-Human

ETFb Antibody

Second Ab: 0.2?g/ml HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

[**IMPORTANT NOTE**]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.