

PAA775Mu01

Polyclonal Antibody to Hexokinase 1 (HK1)

Organism Species: *Mus musculus* (Mouse)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.44mg/mL

UOM: 50µL

Cross Reactivity: Human;Mouse

Applications: WB; IHC; ICC; IP.

[IMMUNOGEN]

Immunogen: Recombinant HK1 (Phe137~Leu270 (Accession # P17710)) expressed in *E.coli*

Accession No.: RPA775Mu01

[APPLICATIONS]

Western blotting: 0.01-2µg/mL;

Immunohistochemistry: 5-20µg/mL;

Immunocytochemistry: 5-20µg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no


obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

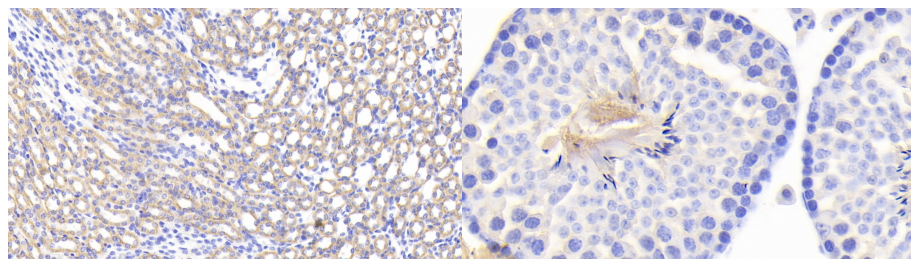


DAB staining on IHC-P; Samples:
Mouse Kidney Tissue; Primary Ab:
20?g/ml Rabbit Anti-Mouse HK1
Antibody Second Ab: 2μg/mL HRP-
Linked Caprine Anti-Rabbit IgG
Polyclonal Antibody (Catalog:
SAA544Rb19)


DAB staining on IHC-P;
Samples: Mouse Testis Tissue;
Primary Ab: 20?g/ml Rabbit Anti-Mouse
HK1 Antibody
Second Ab: 2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

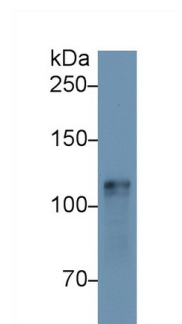


Western Blot; Sample: Mouse Heart
lysate;
Primary Ab: 1μg/ml Rabbit Anti-Mouse
HK1 Antibody
Second Ab: 0.2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

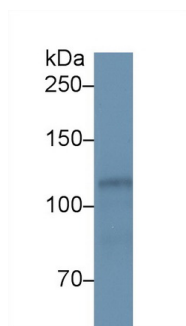


Western Blot; Sample: Mouse Lung
lysate;
Primary Ab: 1μg/ml Rabbit Anti-Mouse
HK1 Antibody
Second Ab: 0.2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

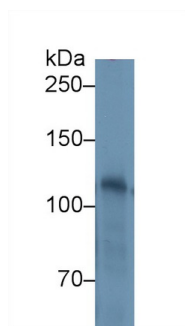

Western Blot; Sample: Mouse Testis
lysate;
Primary Ab: 1μg/ml Rabbit Anti-Mouse
HK1 Antibody
Second Ab: 0.2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

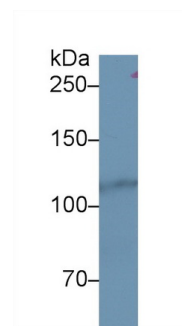

Western Blot; Sample: Human Jurkat
cell lysate;
Primary Ab: 1μg/ml Rabbit Anti-Mouse
HK1 Antibody
Second Ab: 0.2μg/mL HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody
(Catalog: SAA544Rb19)

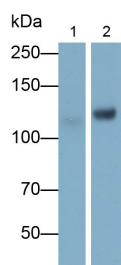

Western Blot; Samples: Lane1: Mouse
Heart lysate; Lane2: Mouse Cerebrum
lysate;

Primary Ab: 0.1 μ g/ml Rabbit Anti-
Mouse HK1 Antibody

Second Ab: 0.2 μ g/ml HRP-Linked
Caprine Anti-Rabbit IgG Polyclonal
Antibody

(Catalog: SAA544Rb19)

[**IMPORTANT NOTE**]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.