

PAA371Ca01

Polyclonal Antibody to Interferon Gamma Induced Protein 10kDa (IP10)

Organism Species: *Canis familiaris*; Canine (Dog)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 200μL

Cross Reactivity: Mouse

Applications: WB

[IMMUNOGEN]

Immunogen: Recombinant IP10 (Ile22~Pro98) expressed in *E.coli*

Accession No.: RPA371Ca01

[APPLICATIONS]

Western blotting: 0.01-2μg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]


Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

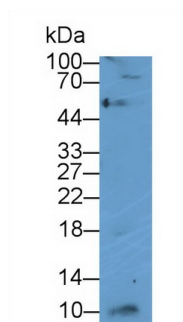

Western Blot; Samples: Lane1: Canine Western Blot; Sample: Mouse YAC1

Ovary lysate; Lane2: YAC1 cell lysate; cell lysate;

Primary Ab: 0.1 μ g/ml Rabbit Anti- Primary Ab: 1 μ g/ml Rabbit Anti-Canine
Canine IP10 Antibody Second Ab: IP10 Antibody

0.2 μ g/ml HRP-Linked Caprine Anti- Second Ab: 0.2 μ g/mL HRP-Linked
Rabbit IgG Polyclonal Antibody Caprine Anti-Rabbit IgG Polyclonal
(Catalog: SAA544Rb19) Antibody

(Catalog: SAA544Rb19)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.