

PAC614Ra01

Polyclonal Antibody to Macrophage Receptor With Collagenous Structure (MARCO)

Organism Species: *Rattus norvegicus* (Rat)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 200μL

Cross Reactivity: Human;Mouse;Porcine

Applications: WB; IHC

[IMMUNOGEN]

Immunogen: Recombinant MARCO (Gly367~Arg519) expressed in *E.coli*

Accession No.: RPC614Ra01

[APPLICATIONS]

Western blotting: 0.01-2μg/mL;

Immunohistochemistry: 5-20μg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the

expiration date under appropriate storage condition.

[IDENTIFICATION]

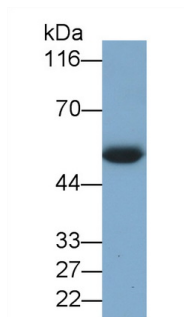

Western Blot; Sample: Rat Thymus lysate


Primary Ab: 0.3µg/ml Rabbit Anti-Rat MARCO Antibody

Second Ab: 0.2µg/mL HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

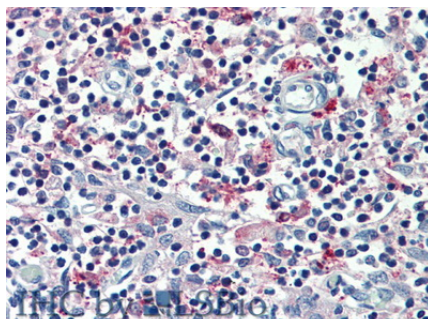

Vector Red staining on IHC-P;


Samples: Human Thymus Tissue;

Primary Ab: 20µg/ml Rabbit Anti-

Human MARCO Antibody Second Ab:

2µg/mL HRP-Linked Caprine Anti-

Rabbit IgG Polyclonal Antibody

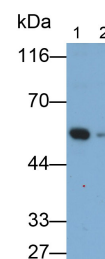

Western Blot; Samples: Lane1: Porcine

Lung lysate; Lane2: Porcine Lymph

node lysate;

Primary Ab: 0.5µg/ml Rabbit Anti-Rat


MARCO Antibody

Second Ab: 0.2?g/ml HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

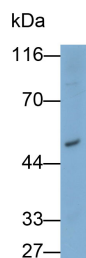

Western Blot; Sample: RAW264.7 cell

lysate

Primary Ab: 0.5µg/ml Rabbit Anti-Rat

MARCO Antibody

Second Ab: 0.2?g/ml HRP-Linked

Caprine Anti-Rabbit IgG Polyclonal

Antibody

(Catalog: SAA544Rb19)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.