

RPX296Hu01 50µg

Recombinant Interleukin 27 (IL27)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Phe29~Pro243+GGGGS*3+Arg21~Lys229

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% Sarcosyl, 5% Trehalose.

Original Concentration: 350µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.0

Predicted Molecular Mass: 52.6kDa

Accurate Molecular Mass: 46kDa as determined by SDS-PAGE reducing conditions.

Phenomenon explanation:

The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affects the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

			FP	RPPGRPQLSL	QELRREFTVS
LHLARKLLSE	VRGQAHFAE	SHLPGVNLYL	LPLGEQLPDV	SLTFQAWRRL	
SDPERLCFIS	TTLQPFHALL	GGLGTQGRWT	NMERMQLWAM	RLDLRDLQRH	
LRFQVLAAGF	NLPEEEEEEE	EEEEEEKGL	LPGALGSALQ	GPAQVSWPQL	
LSTYRLLHSL	ELVLSRAVRE	LLLLSKAGHS	VWPLGFPTLS	PQP	
GGGGSGGGGSGGGGS					
		RKGPPAALT	PRVQCASRY	PIAVDCSWTL	
PPAPNSTSPV	SFIATYRLGM	AARGHSWPCL	QQTPTSTSTCT	ITDVQLFSMA	
PYVLNVTAVH	PWGSSSSFVP	FITEHIIKPD	PPEGVRLSPL	AERQLQVQWE	
PPGSWPFPEI	FSLKYWIRYK	RQGAARFHRV	GPIEATSFIL	RAVRPRARYY	
VQVAAQDLTD	YGELSDWSLP	ATATMSLGK			

[IDENTIFICATION]

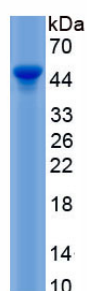


Figure. SDS-PAGE

[IMPORTANT NOTE]



The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.