

RPB166Hu03 200µg

Recombinant VGF Nerve Growth Factor Inducible (VGF)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ala23~Pro615

Tags: N-terminal His Tag

Subcellular Location: Secreted, Cytoplasm

Purity: > 80%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% Sarcosyl, 5% Trehalose.

Original Concentration: 60µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.4

Predicted Molecular Mass: 68.7kDa

Accurate Molecular Mass: 80kDa as determined by SDS-PAGE reducing conditions.

Phenomenon explanation:

The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affects the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

APPGRPEAQPPPLSSEHKPEVAGDAVPGPKDGSAPVVRGARNSEPQDEGELF
 QGVDPRALAAVLLQALDRPASPPAPSGSQQGPEEEAAEALLTETVRSQTHSLP
 APESPEPAAPPRPQTPENGPEASDPSEEELEALASLLQELRDFSPSSAKRQQETA
 AAETETRTHTLTRVNLESPGPERVWRASWGEFQARVPERAPLPPPAPSQFQA
 RMPDSGPLPETHKFGEGVSSPKTHLGEALAPLSKAYQGVAAPFPKARRPESAL
 LGGSEAGERLLQQGLAQVEAGRRQAEATRQAAAQEERLADLASDLLLQYLLQ
 GGARQRGLGGRGLQEAAEERESAREEEEEAEQERRGGEERVGEEDEEAAEAEA
 EAAEAERARQNALLFAEEEDGEAGAEDKRSQEETPGHRRKEAEGTEEGGEEED
 DEEMDPQTIDSLIELSTKLHLPADDVVSIIEEVEEKRRKKNAPPEPVPPPRAAP
 APTHVRSPQPPPPAPAPARDELPDWNEVLPPWDREEDVYPPPGPYHPFPNYI
 RPRTLQPPSALRRRHYYHALPPSRHYPGREAQARRAQEEAAEERRLQEQEEL
 ENYIEHVLLRRP

[IDENTIFICATION]

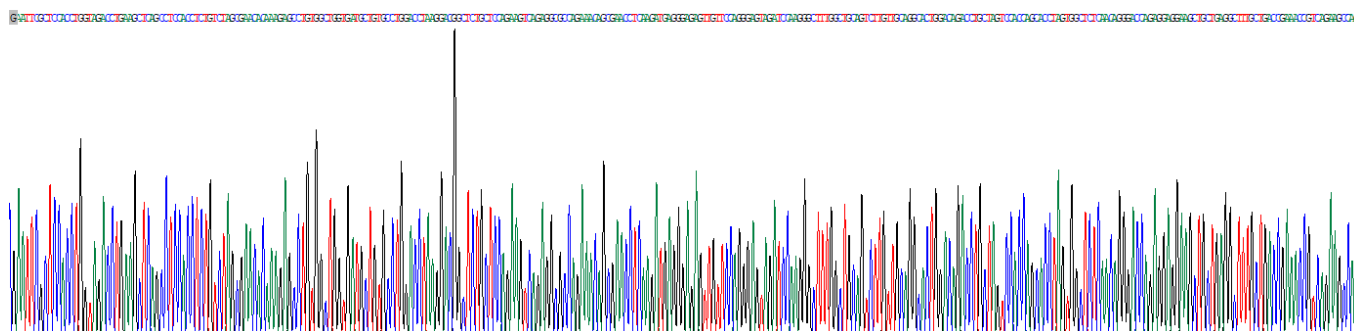


Figure . Gene Sequencing (extract)

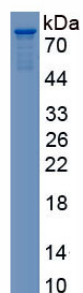


Figure. SDS-PAGE

[**IMPORTANT NOTE**]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.