

**PAE237Mu01****Polyclonal Antibody to Hyaluronoglucosaminidase 2 (HYAL2)****Organism Species: Mus musculus (Mouse)*****Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

**[ PRODUCT INFORMATION ]****Immunogen:** HYAL2, Mouse**Clonality:** Polyclonal**Host:** Rabbit**Immunoglobulin Type:** IgG**Purification:** Affinity Chromatography.**Applications:** WB, ICC, IHC-P, IHC-F, ELISA**Concentration:** 200µg/mL**UOM:** 100µg**[ IMMUNOGEN INFORMATION ]****Immunogen:** Recombinant HYAL2 (Leu23~Gly447) expressed in *E.coli*.**Accession No.:** RPE237Mu01**Sequence:** The target protein is fused with N-terminal His-Tag and its sequence is listed below.

MGHHHHHHHSGSEF-LKPTAPPI FTGRPFVVAW NVPTQECAPR HKVPLDLRAF  
DVKATPNEG F FNQNTTFYY DRLGLYPRFD AAGTSVHGGV PQNGSLCAHL PMLKESVERY  
IQTQEPGGLA VIDWEEWRPV WVRNWQEKDV YRQSSRQLVA SRHPDWPSDR  
VMKQAQYEF E FAARQFMLNT LRYVKAVRPQ HLWGFYLFDP CYNHDYVQNW  
ESYTGRCPDV EVARNDQLAW LWAESTALFP SVYLDETLAS SVHSRNFVSF RVREALRVAH  
THHANHALPV YVFTPTYR GLTGLSQVDL ISTIGESAAL GSAGVIFWGD SEDASSMETC  
QYLKNYLTQL LVPYIVNVSW ATQYCSWTQC HGHGRCVRRN PSANTFLHLN ASSFRLVPGH  
TPSEPQLRPE GQLSEADLNY LQKHFRQCQY LGWGGEQCQR NYKGAAG

**[ ANTIBODY SPECIFICITY ]**

The antibody is a rabbit polyclonal antibody raised against HYAL2. It has been selected for its ability to recognize HYAL2 in immunohistochemical staining and western blotting.

**[ APPLICATIONS ]**

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500


Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200


Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

## [ IMAGES ]

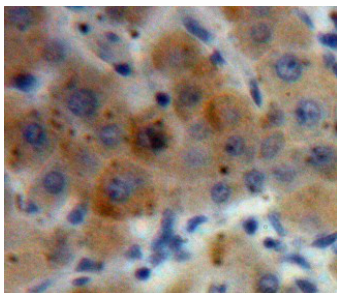

Used in DAB staining on formalin fixed  
paraffin-embedded liver tissue

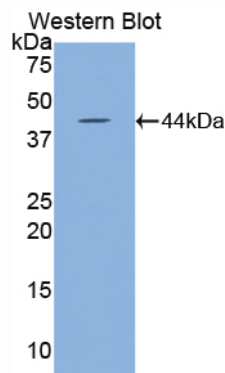

Used in Western Blot, Sample:  
Recombinant HYAL2, Mouse

## [ CONTENTS ]

**Form & Buffer:** Supplied as solution form in PBS, pH 7.4, containing 0.02%  $\text{NaN}_3$ , 50% glycerol.

## [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant HYAL2 (Leu23~Gly447) disposed in loading buffer.

**Usage:** 10  $\mu\text{L}$  per well when 3,3'-Diaminobenzidine (DAB) as the substrate.  
5  $\mu\text{L}$  per well when used in enhanced chemiluminescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris (pH 8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%,  $\text{NaN}_3$  0.02%.

## [ STORAGE ]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.