


A90938Mu01
Polyclonal Antibody to Podocin (PDCN)
Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
 NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

5th Edition (Revised in January, 2013)

[PRODUCT INFORMATION]

Immunogen: PDCN
Clonality: Polyclonal
Host: Rabbit
Species Reactivity: Mouse
Ig type: Rabbit IgG
Purification: Antibodies are purified by target protein affinity chromatography.
Applications: WB, ICC, IHC-P, IHC-F
Form: Liquid
Size: 50μg

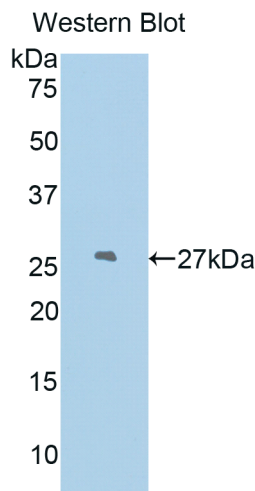

Sample: Recombinant mouse PDCN

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant mouse PDCN (Ile138~Pro371) expressed in *E. coli*.
Molecular Weight: 27.4kDa
USCN accession No.: P90938Mu01
Sequence: The target protein is fused with N-terminal His-Tag and its sequence is listed below.

MGHHHHHHSGSEF-IFR LGHLLPGRAK GPGLFFFLPC LDTYHKVDLR LQTLEIPFHE
 VVTKDMFIME IDAVCYRME NASLLLSSLA HVSKAIQFLV QTTMKRLLAH RSLTEILLER
 KSIAQDVKVA LDAVTCIWGI KVERTEIKDV RLPAGLQHSL AVEAEAQRQA KVRVIAAEGE
 KAASESLRMA AEILSGTPAA VQLRYLHTLQ SLSTEKPATV VLPLPFDMLS LLSSPGNRAQ
 GSINYPSSSK P

Unique product Superb quality Client favorite Nicest service  ISO9001:2008;  ISO13485:2003; 

[**ANTIBODY SPECIFICITY**]

Anti PDCN is a rabbit polyclonal antibody raised against mouse PDCN. This antibody has been selected for its ability to recognize mouse PDCN in immunohistochemical staining and western blotting, non cross-reactive with other members of the family.

[**APPLICATIONS**]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Optimal working dilutions must be determined by end user.

[**CONTENTS**]

Form & Buffer: Supplied as solution form in PBS, pH 7.4, containing 0.02%Na₃N, 50% glycerol.

[**QUALITY CONTROL**]

Content: The quality control contains recombinant mouse PDCN (Ile138~Pro371) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescient (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading buffer: 100mM Tris(PH 8.8), 2%SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, Na₃N 0.02%.

[**STORAGE**]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.