

RPA480Mu01 10μg

Recombinant Myoglobin (MYO)

Organism Species: Mus musculus (Mouse)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Met1~Gly154

Tags: N-terminal His-Tag

Accession: P04247

Host: *E. coli*

Purity: >95%

Endotoxin Level: <1.0EU per 1μg

(determined by the LAL method).

Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 1mM DTT, 5% trehalose, 0.01% sarcosyl and preservative.

Predicted isoelectric point: 7.3


Predicted Molecular Mass: 18.3kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

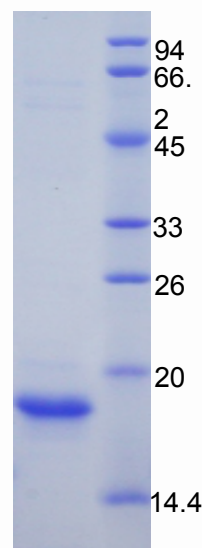

15% SDS-PAGE

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

MGLSDGEWQ LVLNVWGKVE ADLAGHGQEV LIGLFKTHPE TLDKFDKFKN LKSEEDMKGS
EDLKKHGCTV LTALGTILKK KGQHAAEIQP LAQSHATKHK IPVKYLEFIS EIIIEVLKKR
HSGDFGADAQ GAMSKALELF RNDIAAKYKE LGFQG

[REFERENCES]

1. Blanchetot A., *et al.* (1986) Eur. J. Biochem. 159:469-474.
2. The MGC Project Team. (2004) Genome Res. 14:2121-2127.
3. Harris D.E., *et al.* (1985) FEBS Lett. 190:288-292.
4. Wunderlich C., *et al.* (2003) Circ. Res. 92:1352-1358.