

# PAC146Hu01 Polyclonal Antibody to Collagen Type VIII Alpha 1 (COL8a1) Organism Species: Homo sapiens (Human) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

# [ PRODUCT INFORMATION ]

Immunogen: COL8a1, Human

Clonality: Polyclonal

Host: Rabbit

Immunoglobulin Type: IgG


**Purification:** Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

**UOM**: 100μg

#### 9th Edition (Revised in Jul, 2013)

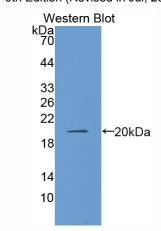

Sample: Recombinant COL8a1, Human

# [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant COL8a1 (Pro575~Pro743) expressed in E.coli.

Accession No.: RPC146Hu01

**Sequence:** The target protein is fused with N-terminal His-Tag and its sequence

is listed below.

MGHHHHHHSGSEF-PPTPPP QGEYLPDMGL GIDGVKPPHA YGAKKGKNGG PAYEMPAFTA ELTAPFPPVG APVKFNKLLY NGRQNYNPQT GIFTCEVPGV YYFAYHVHCKGGNVWVALFK NNEPVMYTYD EYKKGFLDQA SGSAVLLLRPGDRVFLQMPS EQAAGLYAGQ YVHSSFSGYL LYP



#### [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against COL8a1. It has been selected for its ability to recognize COL8a1 in immunohistochemical staining and western blotting.

# [APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

# [CONTENTS]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

# [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant COL8a1 (Pro575~Pro743) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.

# [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.