

PAA097Mu01

Polyclonal Antibody to Matrix Metalloproteinase 1 (MMP1)

Organism Species: *Mus musculus* (Mouse)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 20µL

Cross Reactivity: Human;Porcine;Caprine;Ovine

Applications: WB; IHC

[IMMUNOGEN]

Immunogen: Recombinant MMP1 (Ala98~Pro275) expressed in *E.coli*

Accession No.: RPA097Mu01

[APPLICATIONS]

Western blotting: 0.01-2µg/mL;

Immunohistochemistry: 5-20µg/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the


expiration date under appropriate storage condition.

[IDENTIFICATION]

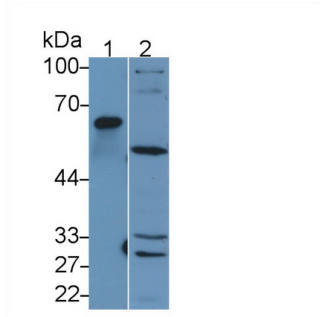

Western Blot; Sample: Lane1: Porcine Cartilage lysate; Lane2: Human A375 cell lysate

Primary Ab: 1µg/mL Rabbit Anti-Mouse MMP1 Ab
Second Ab: 0.2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody
(Catalog: SAA544Rb19)

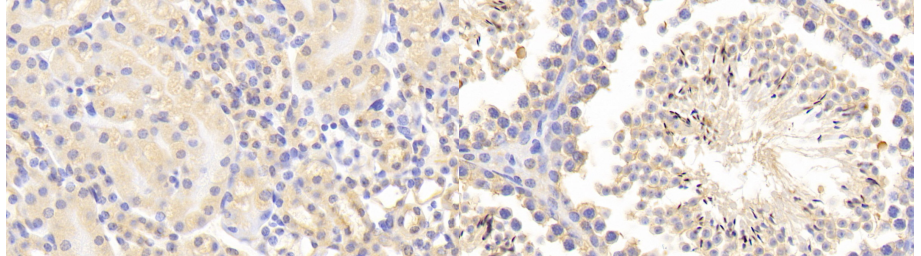


DAB staining on IHC-P; Samples:

Mouse Kidney Tissue; Primary Ab: 20µg/ml Rabbit Anti-Mouse MMP1 Antibody Second Ab: 2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody (Catalog: SAA544Rb19)

DAB staining on IHC-P;


Samples: Mouse Testis Tissue; Primary Ab: 20µg/ml Rabbit Anti-Mouse MMP1 Antibody
Second Ab: 2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody
(Catalog: SAA544Rb19)

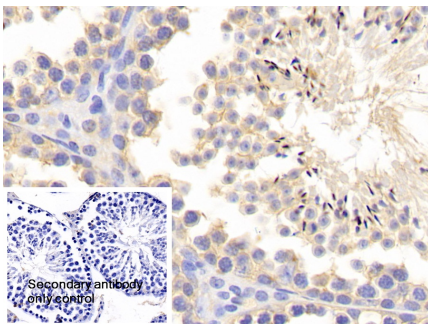

DAB staining on IHC-P;

Sample: Mouse Testis Tissue

Primary Ab: 20µg/ml Rabbit Anti-Mouse MMP1 Antibody
Control: Used PBS instead of primary antibody
Second Ab: 2µg/ml HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody
(Catalog: SAA544Rb19)

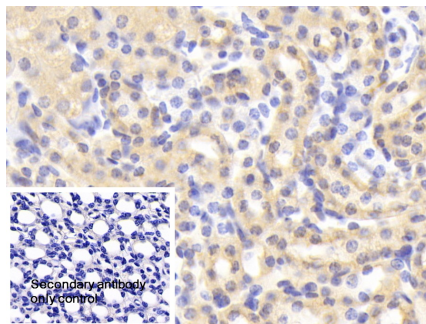

DAB staining on IHC-P;

Sample: Mouse Kidney Tissue

Primary Ab: 20µg/ml Rabbit Anti-Mouse MMP1 Antibody
Control: Used PBS instead of primary antibody
Second Ab: 2µg/ml HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody
(Catalog: SAA544Rb19)



Western Blot; Sample: Caprine Skin lysate

Primary Ab: 0.5µg/ml Rabbit Anti-Mouse MMP1 Antibody
Second Ab: 0.2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody
(Catalog: SAA544Rb19)

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.