

RPA013Hu03 3mg

Recombinant Bone Morphogenetic Protein 2 (BMP2)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Gln283~Arg396

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose .

Original Concentration: 400µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 14.2kDa

Accurate Molecular Mass: 16kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

QAKHKQRK RLKSSCKRHP
 LYVDFSDVGW NDWIVAPPGY HAFYCHGECF FPLADHLNST NHAIVQTLVN
 SVNSKIPKAC CVPTLSAIS MLYLDENEKV VLKNYQDMVV EGCGR

[IDENTIFICATION]

TCAAGCCAAACAAACAGCGGAACCCCTTAAGTCGAGCTGTAAAGAGACACCCCTTTGACGTGGACTTCACTGACGTGGGTGGATGACTGGATTGTGCTCCCGGGGTATCACGCTTTTACTGCGACGAGATGCCCTTTCTGCTGCTGATCATCTGAACCTCACTAATCATGCCATTGTTCAAGACGTTGGTCACTCTGTT
 Q A K H K Q R K R L K S S C K R H P L Y V D F S D V G W H D W I V A P P G Y H A F Y C H G E C F F P L A D H L N S T N H A I V Q T L V N S V I

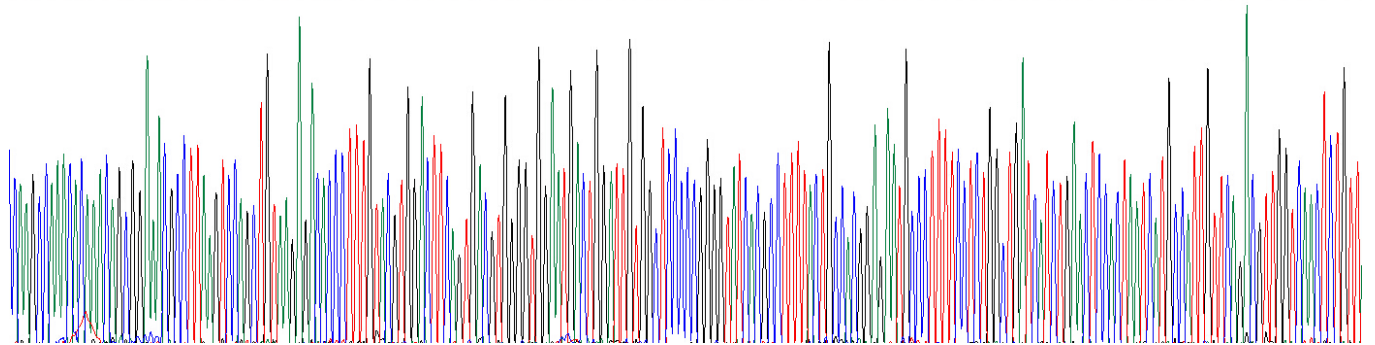


Figure . Gene Sequencing (extract)

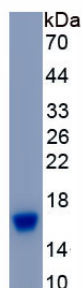


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.