

RPF233Hu01 100µg

Recombinant Brain Specific Angiogenesis Inhibitor 3 (BAI3)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Cys30~Pro508

Tags: Two N-terminal Tags, His-tag and T7-tag

Accession: O60242

Host: *E. coli*

Subcellular Location: Multi-pass membrane protein, Cell membrane.

Purity: >90%

Endotoxin Level: <1.0EU per 1µg
(determined by the LAL method).


Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 5% trehalose, 0.01% sarcosyl.

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 57.9kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

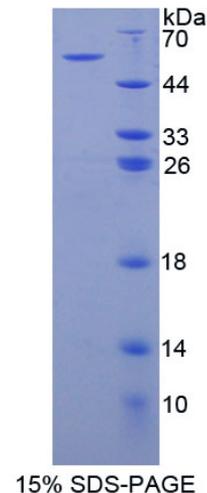

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

C STLVKGVYIG SYSVSEMFPK NFTNCTWTLE NPDPTKYSIY LKFSKDLSC SNFSLLAYQF
DHFSHEKIKD LLRKNHNSIMQ LCNSKNAFVF LQYDKNFIQI RRVFPTNFPG LQKKGEEDQK
SFFEFLVLNK VSPSQFGCHV LCTWLESCLK SENGRTESCG IMYTKCTCPQ HLGEWGIDDQ
SLILLNNVVL PLNEQTEGCL TQELQTTQVC NLTREAKRPP KEEFGMMGDH TIKSQRPRSV
HEKRVPQEQA DAAKFMAQTG ESGVEEWSQW STCSVTCGQG SQVRTRTCVS
PYGTHCSGPL RESRVCNNTA LCPVHGVWEE WSPWSLCSFT CGRGQRTRTR
SCTPPQYGGR PCEGPETHHK PCNIALCPVD GQWQEWSSWS QCSVTCSNGT
QQRSRQCTAA AHGGSECRGP WAESRECYNP ECTANGQWNQ WGHWSGCSSK
CDGGWERRIR TCQGAUITGQ QCEGTGEEVR RCNEQRCP