

RPA853Hu03 100ug

Recombinant Caspase 8 (CASP8)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Leu7~Asp201

Tags: N-terminal His Tag

Subcellular Location: Chromosome

Purity: > 97%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.8

Predicted Molecular Mass: 26.6kDa

Accurate Molecular Mass: 28kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

LYDI GEQLDSEDLA SLKFLSLDYI PQRKQEPIKD ALMLFQRLQE KRMLEESNLS FLKELLFRIN RLDLLITYLN
 TRKEEMEREL QTPGRAQISA YRVMLYQISE EVSRSELRSF KFLLQEEISK CKLDDDMNLL DIFIEMEKRV ILGEGKLDIL
 KRVCAQINKS LLKIINDYEE FSKGEELCGV MTISDSPREQ D

[IDENTIFICATION]

DCCTTATGATATGGGGNCAAC TG GACAGTGAAGACTTGGCTCCCTCAGTTCCTGAGCCTGGACTACATTCGCAAGGAAGCAAGAACCAACAAGSATGCCITGATGTTATTCAGAGACTCCAGGAAAGAGAAATGTGGAGGAAGCAATCTGTCTCTCCGAGGAGCTGCTCTCCGAAITAAAGACTGGATTTGCTGATTAC
 L Y D I G E Q L D S E D L A S L K F L S L D Y I P Q R K Q E P I K D A L M L F Q R L Q E K R M L E E S N L S F L K E L L F R I N R L D L L I T


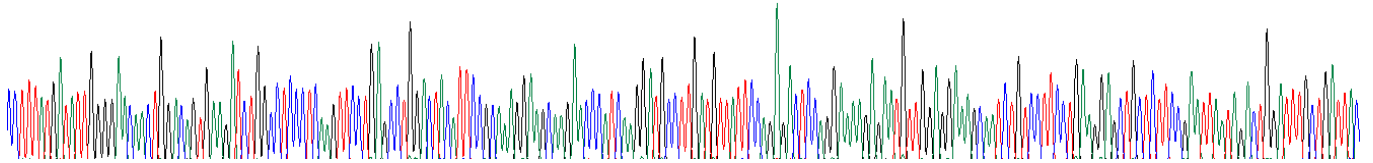


Figure. Gene Sequencing (Extract)

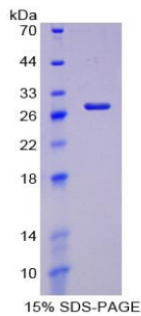

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.