

RPB882Hu01 50µg

Recombinant Fibroblast Growth Factor 10 (FGF10)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Gln38~Ser208

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% skl, 5%Trehalose.

Original Concentration: 450µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 10.0

Predicted Molecular Mass: 23.1kDa

Accurate Molecular Mass: 23kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

QAL QDMVSPEAT
NSSSSSFSSP SSAGRHVRSY NHLQGDVRWR KLFSFTKYFL KIEKNGKVSG
TKKENCYPSI LEITSVEIGV VAVKAINSNY YLAMNKKGKL YGSKEFNNDK
KLKERIEENG YNTYASFNWQ HNGRQMYVAL NGKGAPRRGQ KTRRKNTSAH
FLPMVVHS

[IDENTIFICATION]

CGATTATGAGCCCTGGTGGCATGATGGAGGACGACCTCTCTCCCTCCTTCCTCTCTGAGCGGGGAGGATGGCGGCTGATCCCTGAGGATGCTGGGAGTATCTCTTACCGATCTTCTGAGTGGAGAGCGGAGGCTGGGACCGAGGAGACTGGCTGACATCTGGGATGATCTGAT

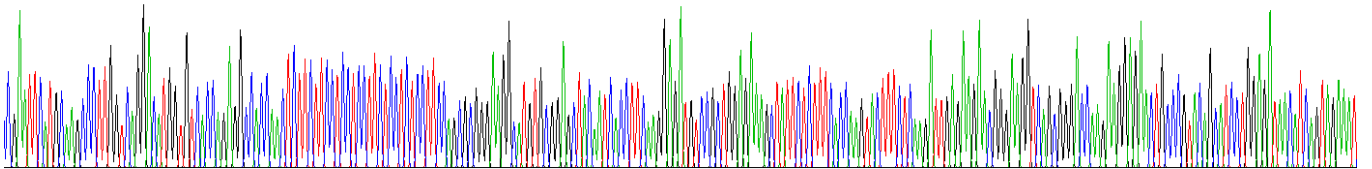


Figure . Gene Sequencing (extract)

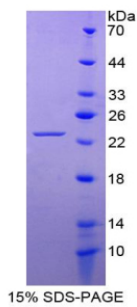


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.