

RPA212Hu01 10µg

Recombinant Matrix Metalloproteinase 23B (MMP23B)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Tyr79~Gly390

Tags: N-terminal His Tag

Subcellular Location: n/a

Purity: > 97%

Traits: Freeze-dried powder

Buffer formulation: 100mMNaHCO₃, 500mMNaCl, pH8.3, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.6

Predicted Molecular Mass: 39.8kDa

Accurate Molecular Mass: 40kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 100mM NaHCO₃, 500mM NaCl (pH8.3) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

YT LTPARLRWDH FNLTyrILSF
PRNLLSPRET RRALAAAFRM WSDVSPFSFR EVAPEQPSDL RIGFYPIHNT
DCLVSALHHC FDGPTGELAH AFFPPHGGIH FDDSEYWVLG PTRYSWKKGV
WLTDLVHVAA HEIGHALGLM HSQHGRALMH LNATLRGWKA LSQDELWGLH
RLYGCLDRLF VCASWARRGF CDARRRLMKR LCPSSCDFCY EFPFPTVATT
PPPPRTKTRL VPEGRNVTFR CGQKILHKKG KVYWKDQEP LEFSYPGYLA
LGEAHLIIA NAVNEGTYTC VVRRQQRVLT TYSWRVVRVG

[IDENTIFICATION]

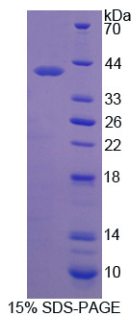

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.