

**RPP155Ra01 10µg**

**Recombinant Wingless Type MMTV Integration Site Family, Member 3A (WNT3A)**

**Organism Species: *Rattus norvegicus* (Rat)**

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Met1~Phe297

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted, Extracellular matrix

**Purity:** > 95%

**Traits:** Freeze-dried powder

**Buffer formulation:** 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% skl, 5%Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 7.8

**Predicted Molecular Mass:** 41.0kDa

**Accurate Molecular Mass:** 41kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

MYNLCVVVVV VAEVLTQQLN DRALTSRVLL WSLAVGPQYS SLSTQPILCA  
SIPGLVPKQL RFCRNYVEIM PSVAEGVKAG IQECQHQFRG RRWNCTTVSN  
SLAIFGPVLD KATRESAFVH AIASAGVAFA VTRSCAEGSA AICGSSRLQ  
GSPGEGWKWG GCSEDIIEFGG MVSREFADAR ENRPDARSAM NRHNNEAGRQ  
AIASHMHLKC KCHGLSGSCE VKTCWWSQPD FRTIGDFLKD KYDSASEMVV  
EKHRESRGWV ETLRPRYTYF KVPTERDLVY YEASPNFCEP NPETGSF

**[ IDENTIFICATION ]**

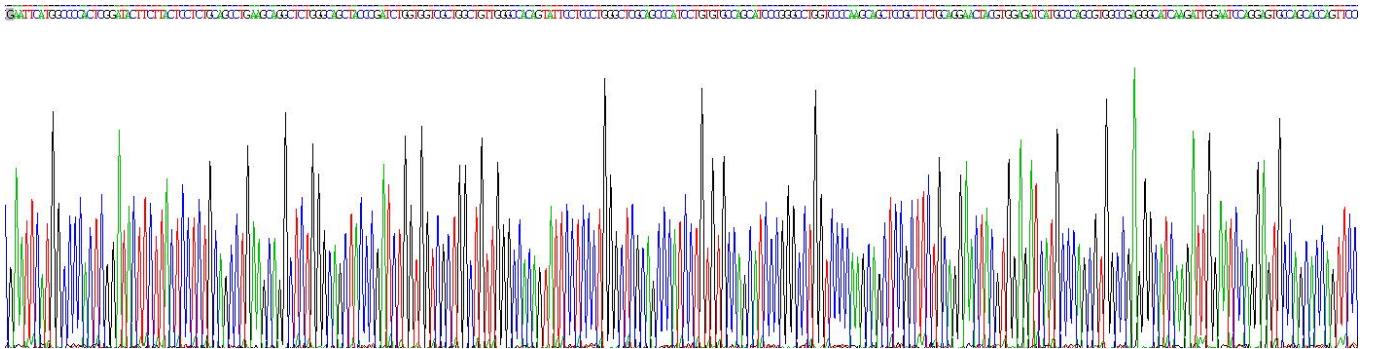

**Figure . Gene Sequencing (extract)**

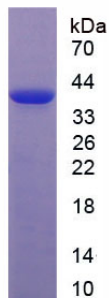


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.