

**RPD231Hu01 10µg**  
**Recombinant Troponin T Type 1, Slow Skeletal (TNNT1)**  
**Organism Species: Homo sapiens (Human)**  
***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

## **[ PROPERTIES ]**

**Residues:** Met1~Gly259 (Accession # P13805),  
with N-terminal His-Tag.

**Host:** *E. coli*

**Purity:** >95%

**Endotoxin Level:** <1.0EU per 1µg  
(determined by the LAL method).

**Formulation:** Supplied as lyophilized form in PBS,  
pH7.4, containing 5% sucrose, 0.01% sarcosyl.

**Predicted isoelectric point:** 5.7

**Predicted Molecular Mass:** 32.3kDa


**Accurate Molecular Mass:** 40kDa as determined by SDS-PAGE reducing conditions.

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

**Note:** The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affects the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

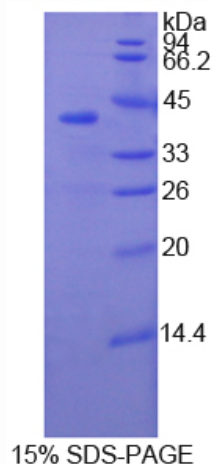

## [ USAGE ]

Reconstitute in sterile PBS, pH7.2-pH7.4.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

## [ SEQUENCES ]

The target protein is fused with N-terminal His-Tag, its sequence is listed below.  
MGHHHHHSGSEF- MSDTEEQEYE EEQPEEEEAE EEEEEAPEEPE PVAEPEEERP  
KPSRPVVPPL IPPKIPEGGER VDFDDIHRKR MEKDLLELQT LIDVHFEQRK KEEELVALK  
ERIERRRSER AEQQRFRTEK ERERQAKLAE EKMRKEEEEA KKRAEDDAKK KKVLSNMGH  
FGGYLVKAEQ KRGKRQTGRE MKVRILSERK KPLDIDYMGE EQLREKAQEL SDWIHQLESE  
KFDLMAKLKQ QKYEINVLYN RISHAQKFRK GAGKGRVGG

## [ REFERENCES ]

1. Gahlmann R., *et al.* (1987) J. Biol. Chem. 262:16122-16126.
2. Samson F., *et al.* (1994) Biochem. Biophys. Res. Commun. 199:841-847.
3. Barton P.J.R., *et al.* (1999) Genomics 57:102-109.
4. Novelli G., *et al.* (1993) Cell Biochem. Funct. 11:187-191.